


BiopC[®]
Interventional

full-sCore[®]

**+++ NEW +++
available in 18G**




Where **Medicine**
meets **Engineering**

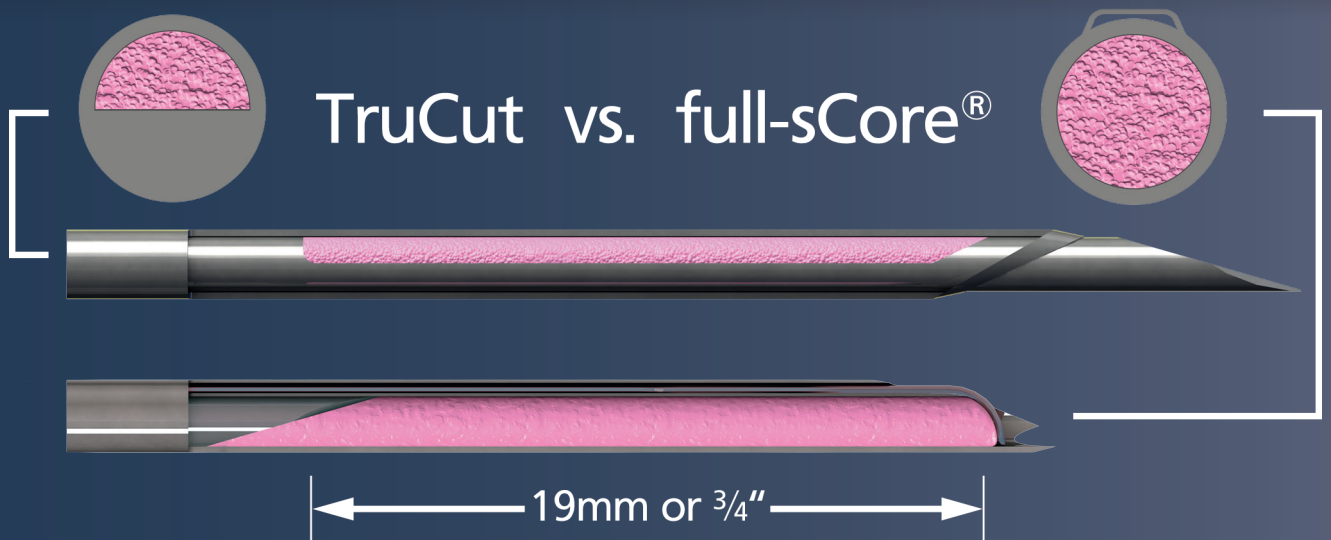


How full-sCore® works


- Advance full-sCore® device to biopsy site

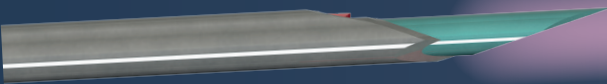


- Gun is triggered. Exterior cannula cuts into tissue.

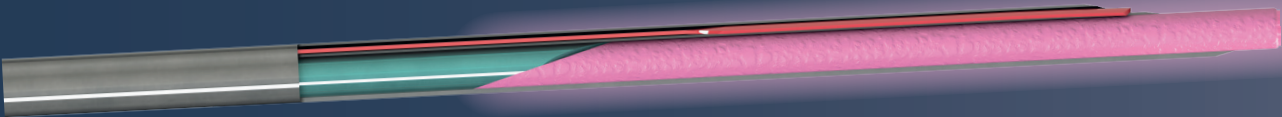

- Instantly, the interior chisel cuts off the tissue sample. Done.

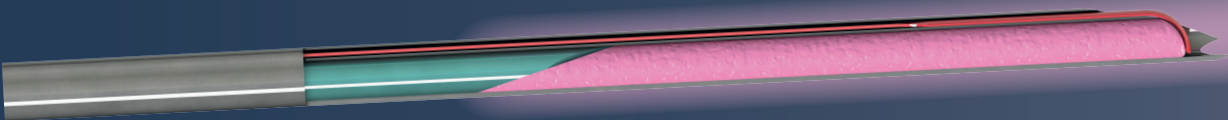

Where Medicine
meets Engineering

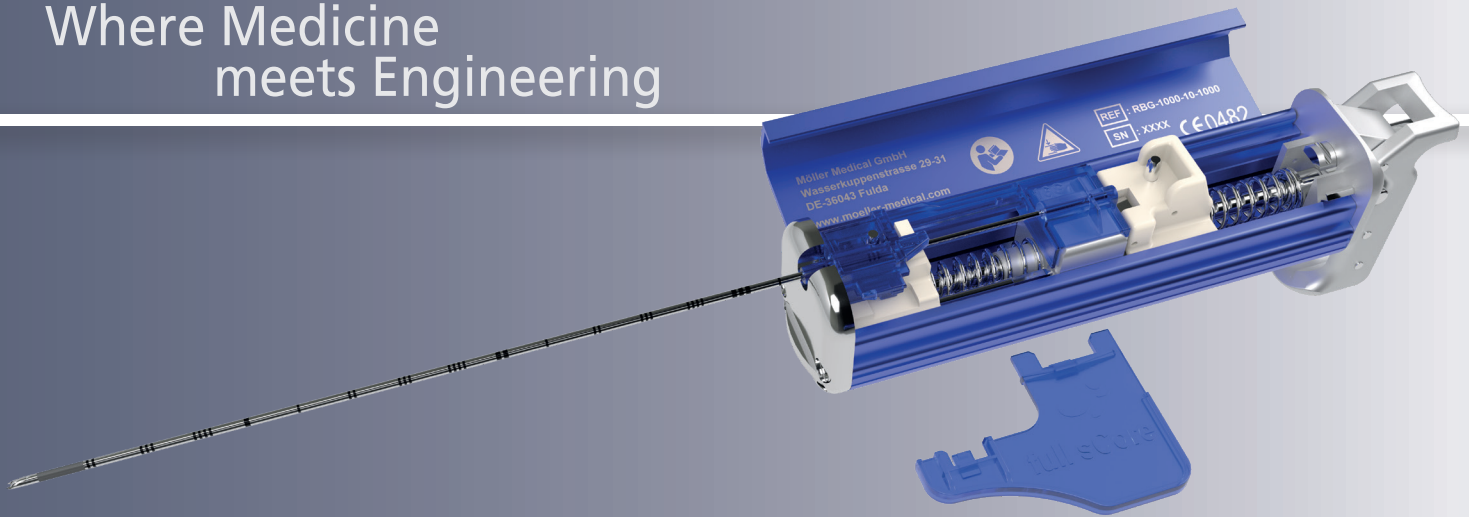

full-sCore® Features


- Möller Medical superfinish provides an optimal surface
- Well-adjusted crown-core cut ensures less traumatic puncture and accurate collection of tissue
- The needle comes with cm markings for a safe and appropriate determination of the depth placement
- An echogenic distal tip ensures superb visibility under ultrasound guided biopsy
- Compatible with COX-2020 coaxial introducer needle for guiding and for achieving multiple sample collection with a single puncture

full-sCore® DNG-2020

Biopsy needle gauge	stitch length	Biopsy needle catalogue ref.	PU
18G	10 cm	DNG-2020-18-0100	10
18G	16 cm	DNG-2020-18-0160	10
18G	20 cm	DNG-2020-18-0200	10
18G	25 cm	DNG-2020-18-0250	10
16G	10 cm	DNG-2020-16-0100	10
16G	16 cm	DNG-2020-16-0160	10
16G	20 cm	DNG-2020-16-0200	10
16G	25 cm	DNG-2020-16-0250	10

COX-2020

stitch length	catalogue ref.	PU
4.9 cm	COX-2020-18-0049	10
10.9 cm	COX-2020-18-0109	10
14.9 cm	COX-2020-18-0149	10
4.9 cm	COX-2020-16-0049	10
10.9 cm	COX-2020-16-0109	10
14.9 cm	COX-2020-16-0149	10

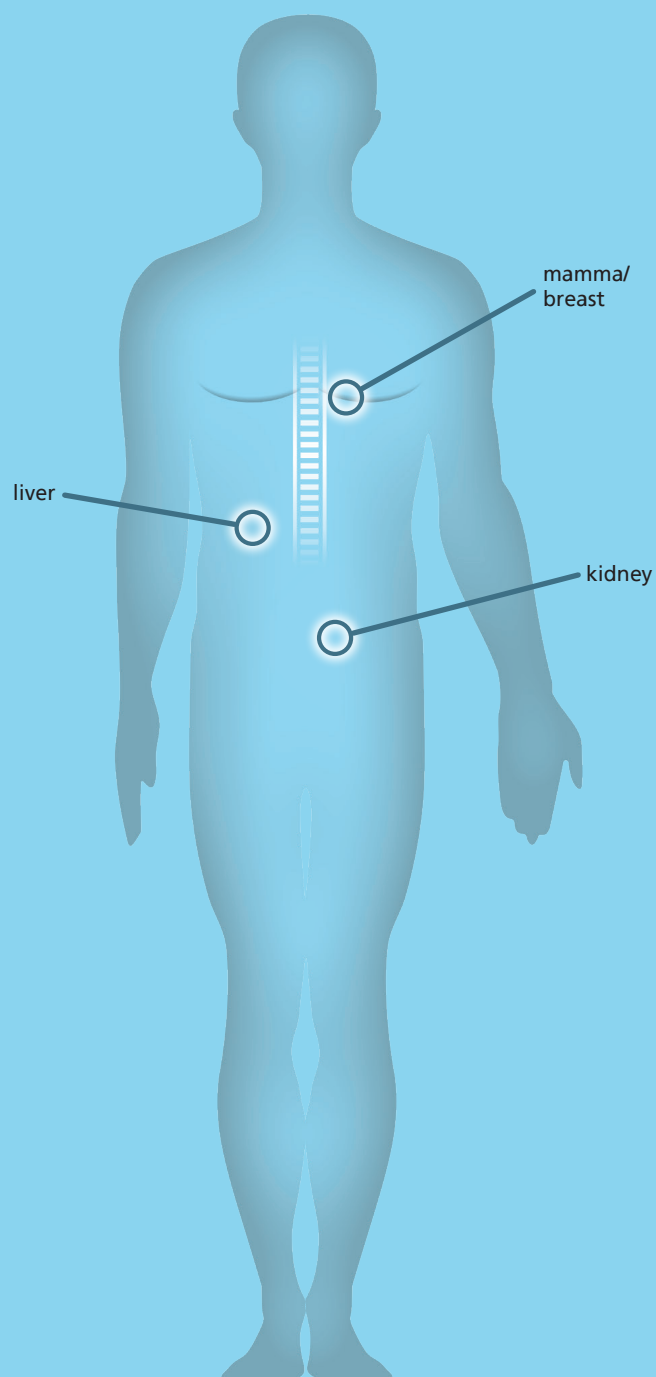

MöllerMedical

Möller Medical GmbH
Wasserkuppenstraße 29-31
D - 36043 Fulda, Germany
Fon +49 661 9 41 95 - 0
Fax +49 661 9 41 95 - 8 50
info@moeller-medical.com
www.moeller-medical.com